


Licking unit

OPR-LKR



- Spout : Sonde
- Sensor : IR photo beam sensor
- Position : XZ manipulator(range of stroke:±15mm)
- Amp box : DC5V driven
/ photo-beam amp circuit
/ 5 V TTL signal output

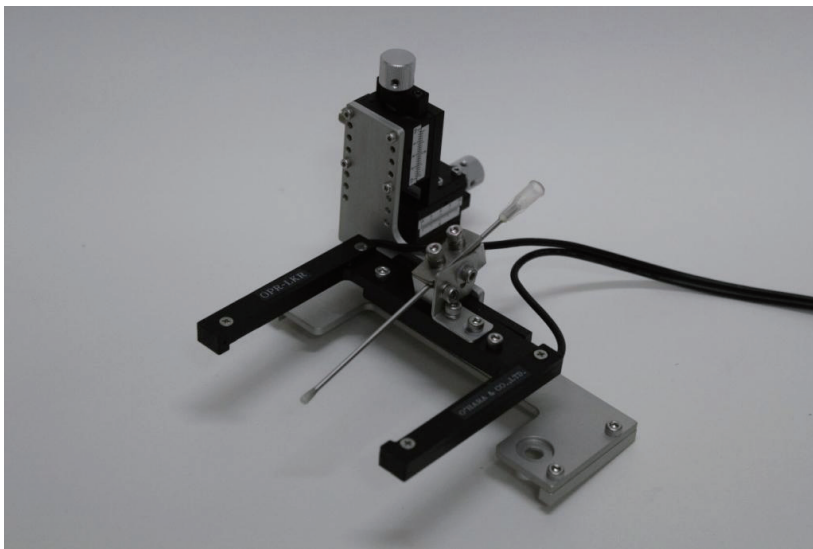

Licking sensor with XYZ manipulator

This is detecting sensor unit that subject animal licks spout. It uses infrared photo beam sensor. It is available to combine with our water supply unit.

★ Please understand that appearance & specification of this product may be changed for improvement without notice.

manufacturer



O' HARA & CO., LTD.

4-28-16 Ekoda, Nakano-ku, Tokyo 165-0022, JAPAN

TEL 81-3-3389-2451

FAX 81-3-3389-2453

<http://www.ohara-time.co.jp>