

Food deprivation tool

ONE
for
10 cages

Automatic food deprivation is available!



Feeder unit

Timer box



This tool is for automatic food deprivation for subject. The subject can approach and eat food only specified period that you set by the timer box.

Using this tool helps to reduce time and effort in performing food deprivation. You can set to most of breeding cage. Parts of the tool that are touched by the subject can all be washed and autoclaved.




Easy
connection
by
only
one cable

Digital type
timer switch,
ON/OFF control
for
12 sets/day
on
1min basis




Food deprivation tool for 10 cages

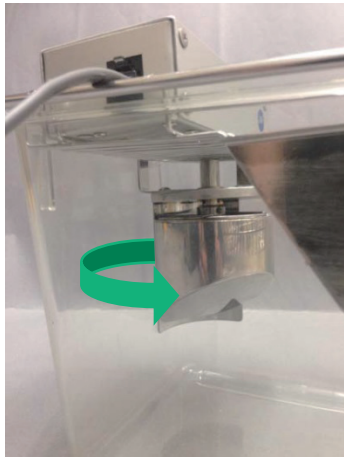

Deprivation

Animal can not approach to foods,
because a cover hides slits of a feedbox.

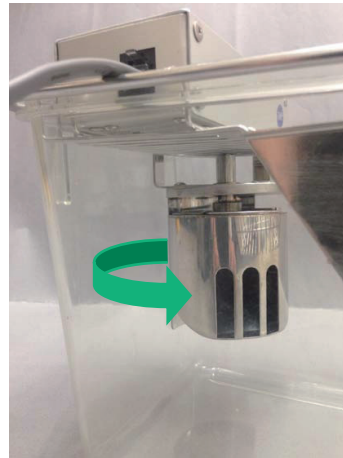

Feed

Animal can approach to foods,
because slits of a feed box face to front.

The motor rotates 180 degree by 2 sec. so it will not hurt or pinch subject



By using
special clasps,
you can use
your own cage!

The feed box is autoclavable.

manufacturer



O' HARA & CO., LTD.

4-28-16 Ekoda, Nakano-ku, Tokyo 165-0022, JAPAN

TEL 81-3-3389-2451

FAX 81-3-3389-2453

<http://www.ohara-time.co.jp>