


Odor Supply Unit

SD-2060

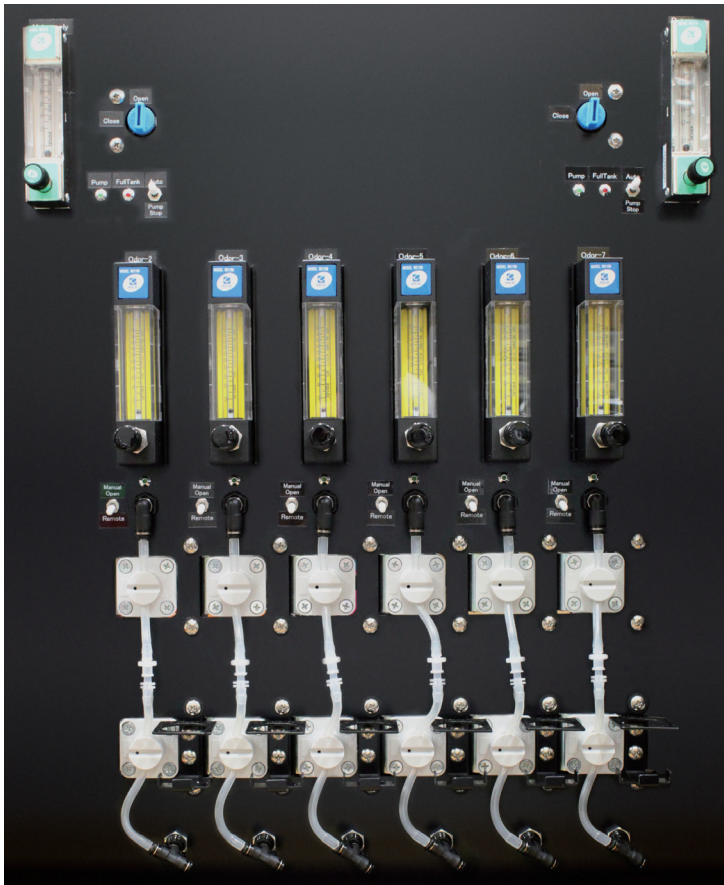

- Odor supply volumes
: 150mL/min
- Main supply and Return pump volumes
: 2000ml/min
- Source of Odor
: Glass vial or Filter mesh

Using this unit, 6 kinds of odors can be presented. Odor type and duration can be controlled programmatically through 5V TTL signal. Air at installation place is filtered and used by this unit. Normally, only the filtered air flows in front of animal's nose. Based on TTL signal, electro magnetic valve opens and odor gets mixed into flowing air. Odor source can be placed in syringe filter or odor bottle. The exhaust air is released to atmosphere flowing through water & an activated carbon filter.




O' HARA & CO., LTD.

★ Please understand that appearance & specification of this product may be changed for improvement without notice.

Odor Supply Unit

SD-2060

System Diagram

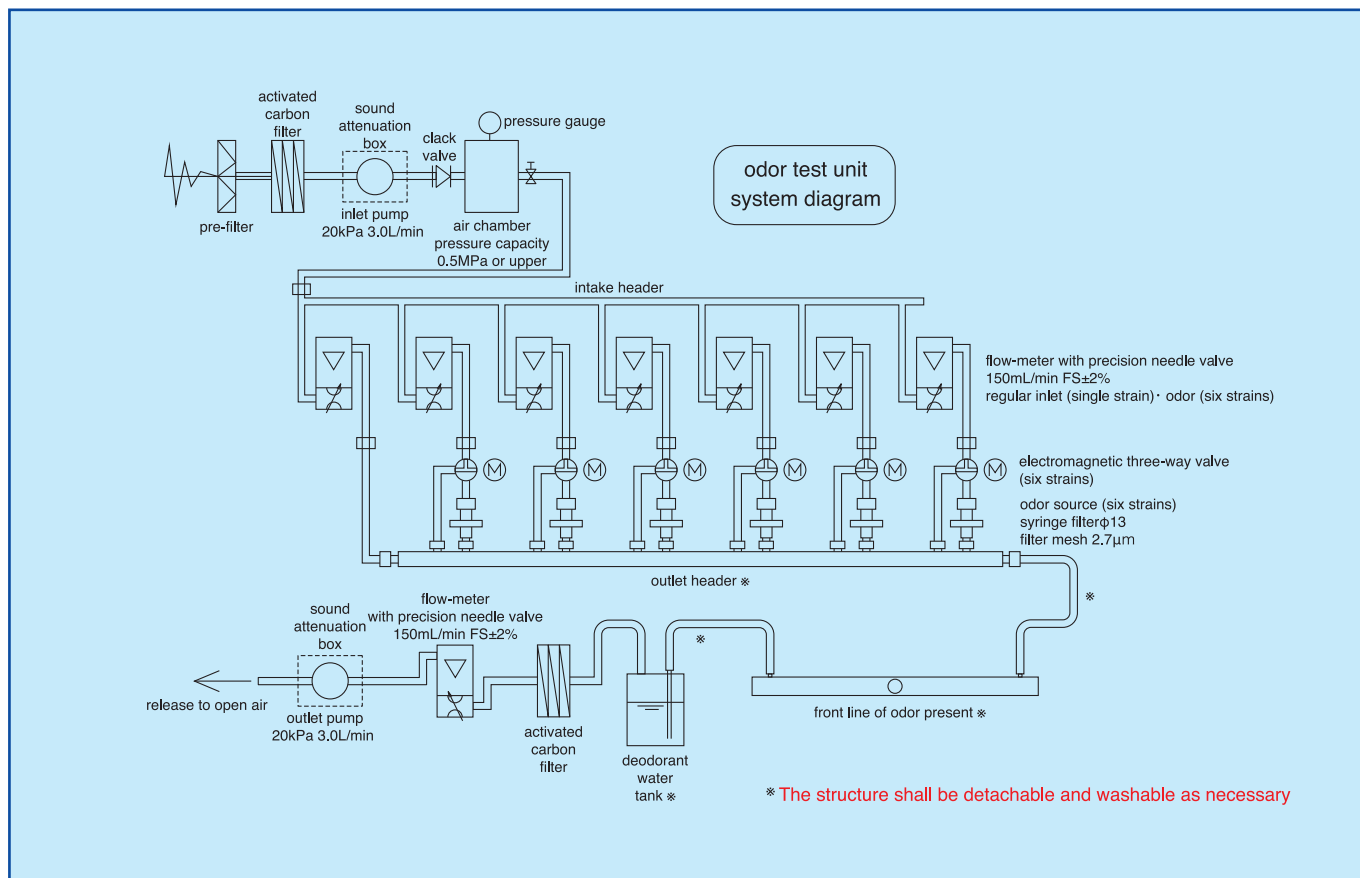

System components

- Volume control panel
- Odor supply unit
- Air supply and vacuum pump unit
- Activated carbon filter unit
- Control interface

manufacturer



O' HARA & CO., LTD.

4-28-16 Ekoda, Nakano-ku, Tokyo 165-0022, JAPAN

TEL 81-3-3389-2451

FAX 81-3-3389-2453

<http://www.ohara-time.co.jp>